


Billerud type pneumatic sample cutter for ECT

Product code [RL-ECT-A](#)

RYCOLAB

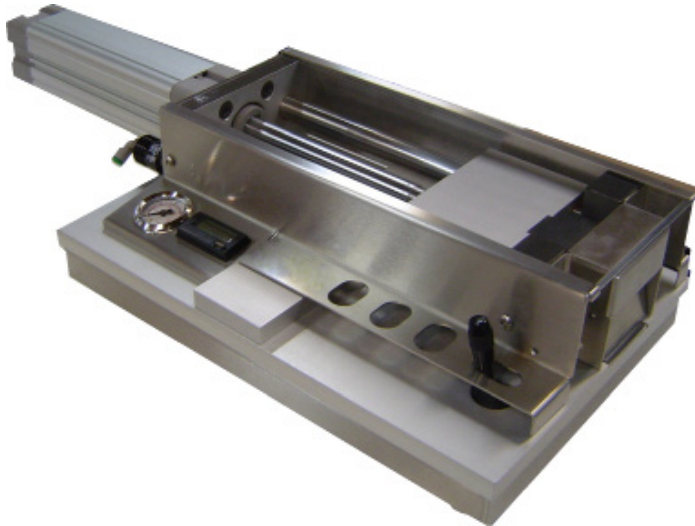

Preparation of corrugated board samples for ECT test.

Applicable standards

ISO 3037, TAPPI T838, T839, etc.

Characteristics

- Pneumatic Cutting that guarantees clean and parallel cuts.
- Digital cuts counter for planning blades substitutions.
- Cutting system with two parallel blades (blades with one edge, to obtain samples without lateral deformations).
- Automatic sample rejection.
- Easy operation equipment.
- Special sample inlet that facilitates the introduction of the sample to be tested.
- Methacrylate protection cover.
- 2 models available : Samples of 25 mm (CB-25) or 2 inches (CB-2i).
- Maximum sample thickness : 15 mm.

Connections:

- Air supply: Max 600 kPa

Weight and dimensions

610 x 250 x 140mm (W x L x H) / 10 kg (only equipment)
720 x 520 x 320mm (W x L x H) / 20 kg (with packing box)

